# Specimen Collection and Submission Guidance for Measles (Rubeola) Serum Specimens

Submitting Serum Specimens for Detection of IgG/IgM Antibodies to Measles

### **Measles Serum Specimen Collection and Storage**

**Report suspected measles cases immediately.** Reporting contacts by county/public health region available <a href="here">here</a>.

Required Specimen: Serum

Required Volume: 200 µL, but more is preferred

#### **Required Storage and Shipping Temperatures:**

- Store and ship cold at 2°C-8°C if specimen will arrive at the Laboratory within 48 hours of collection.
- **Store** and **ship frozen** at -20°C or colder if specimen will arrive at Laboratory **more than 48 hours** after collection.

Ensure specimen collection kits are not expired!

### **Specimen Shipping and Labeling Requirements**

**Ship as:** Category B Biological Substance, UN3373 Specimens must be:

- **Triple packaged** to withstand shock, pressure changes, leaks, and other ordinary handling conditions while in transit.
- Packaged with enough absorbent material in secondary container to soak up the contents of the primary container.
- **Shipped overnight** in insulated containers with frozen cold packs (for cold specimens), or dry ice (for frozen specimens).

**Ensure** containers are securely closed to prevent leaks.

- **Secure** specimen lids by wrapping in paraffin film (e.g., Parafilm).
- Pack cold specimens with multiple ice packs.

Visit DSHS' online **Specimen Shipping and Mailing Guidance** for more details on shipping Category B substances.

Ensure "Saturday Delivery" box is selected for Friday shipments. Do not ship for Sunday or holiday delivery.

Specimens received out of temperature range will be rejected.



A child's cheek showing the characteristic measles rash. (CDC/Tatiana Lanzieri, MD, MPH, 2024)





## **Specimen Collection and Submission Guidance for**

### Measles (Rubeola) Swab Specimens

### **Submitting Measles Swab Specimens for Viral RNA Detection**

### **Measles Swab Specimen Collection and Storage**

**Report suspected measles cases immediately.** Reporting contacts by county/public health region available here.

#### **Required Specimens:**

- Preferred: Throat swabs
- Acceptable: Nasopharyngeal or oral swabs
   Use dacron or rayon swabs only; no cotton tips or wooden shafts.
- Use Universal Transport Media or Viral Transport Media.

### Required Volume: 1 mL-2 mL

### **Required Storage and Shipping Temperature:**

- **Store** and **ship cold** at 2°C—8°C if specimen will arrive at Laboratory **within 48 hours** of collection.
- Store and ship frozen at -70°C if specimen will arrive at Laboratory more than 48 hours after collection.

Ensure specimen collection kits are not expired!

### **Specimen Shipping and Labeling Requirements**

**Ship as:** Category B Biological Substance, UN3373 Specimen must be:


- **Triple packaged** to withstand shock, pressure changes, leaks, and other ordinary handling conditions while in transit.
- Packaged with enough absorbent material in secondary container to soak up the contents of the primary container.
- **Shipped overnight** in insulated containers with enough cold packs to keep them cold or enough dry ice to keep them frozen.
  - Pack cold specimens with multiple frozen ice packs.
- Ensure containers are securely closed to prevent leaks.
  - **Secure** lids by wrapping in paraffin film (e.g., Parafilm).

Visit DSHS' online <u>Specimen Shipping and Mailing Guidance</u> for more details on shipping Category B substances.

Ensure "Saturday Delivery" box is selected for Friday shipments.

Do not ship for Sunday or holiday delivery.

Specimens received out of temperature range will be rejected.

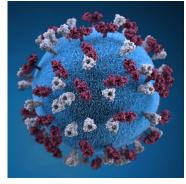

(CDC/ Allison M. Maiuri, MPH, CHES (2016))



# **Specimen Collection and Submission Guidance for**

**Measles Serum Specimens for Antibody Detection** \*\*\*All Facilities Must Have a DSHS Submitter Account to Submit Specimens\*\*\*

### **Label Specimen With Unique Identifiers**

Every specimen must have at least two unique patient identifiers on its label.




Three patient identifiers provided on this label:

- 1. Name
- 2. **Date of Birth**
- 3. **Medical Record Number**

#### Provide Patient Identifiers in Sections 2 and 3 of Form G-2A

Patient identifiers on specimen label and G-2A submission form must match.

Date of Collection must be provided in Section 3.

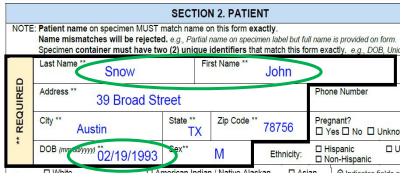

	SECTION 3. SPECIMEN			
NOTE: If the 'Date of Collection' field is not completed, the specimen will be rejected.				
	Date of Collection (mm/dd/yyyy) ** 02/21/2025	Time of Collecti	on ** □ AM □ PM	Col
ED	Unique Identification Number ** e.g., MRN / Alien # / Accession I 06	161858	Comments or Additio e.g., CDC ID, Previous DSH	

## **Identify Reason for Submission in Section 2**

Check "Outbreak" or "Surveillance".

☐ Outbreak Association Θ	
☐ Surveillance Θ	

### **Questions About ...**

Specimen Collection/Suitability: 512-776-7594 or 512-776-7760

**Specimen Shipping:** 512-776-7598 or 1-888-963-7111 ext. 7578 (toll free) **Measles Surveillance:** EAIDUMeasles2025@dshs.texas.gov and Measles (Rubeola) | Texas DSHS

Submitter Accounts, Submission Forms, or Result Reports: 512-776-7578 or LabInfo@dshs.texas.gov



Check "Serum."

"Measles IgG."

One or both tests

may be selected.

**Check** the appropriate Payor. Do not leave empty.

ED	☐ Medicaid (2)  Medicaid/Medicare #:	☐ Medicare (8)
** REQUIRED	□ Submitter (3) □ BIDS (1720) □ BT Grant (1719) □ HIV / STD (1608) □ IDEAS (1610)	☐ Immunizations (1609) ☐ Private Insurance* (4) ☐ TIPP (5144) ☐ Zoonosis (1620) ☐ Other:

**Select Specimen Type in Section 3** 

**Select Test Type in Section 4** 

Check "Measles IgM" and/or



SECTION 6. SEROLO

☐ Brucella, Total Antibody

☐ Hantavirus IgM & IgG ☐ Measles IaM ☐ Measles lαG

☐ Chagas IgG

## **Specimen Collection and Submission Guidance for**

Measles Swab Specimens for Viral RNA Detection

\*\*\*All Facilities Must Have a DSHS Submitter Account to Submit Specimens\*\*\*

### **Label Specimen With Unique Identifiers**

Every specimen must have at least two unique patient identifiers on its label.



- 2 DOB: 02/19/1993
- 3 06161858



- 2. Date of Birth
- 3. Medical Record Number

### Provide Patient Identifiers in Sections 2 and 3 of Form G-2V

**Patient identifiers** on specimen label and G-2V submission form **must match**.

**Date of Collection must** be provided in Section 3.

		SEC	TION 2. PATI	ENT		
NOTE	: Patient name on specimen MUST r Name mismatches will be rejected Specimen container must have tw	d. e.g., Pa	artial name on spe	cimen label but fu		
** REQUIRED	Last Name ** Snow		First Name **	John	>	
	Address ** 39 Broad Str	Broad Street			Phone Number	
	City ** Austin	State *	Zip Code	78756	Pregnant? ☐ Yes ☐ No ☐ Unkn	
	DOB (mm/ad/yyyy) ** 02/19/1993	Sex**	М	Ethnicity:	☐ Hispanic ☐ Non-Hispan	ic U
□ White □ American Indian / Native Alackan □ Asian □ Indicator folds a						
			SECTIO	N 3. SPECI	MEN	
NOT	E: If the 'Date of Collection' fiel	d is not	completed, th	e specimen w	/ill be rejected	
	Date of Collection (mm/dd/yy	yy) **	Time of Colle	ection **	=	Col
	02/21/2025	>			□ AM □ PM	
ΞD	Unique Identification Numbe e.g., MRN / Alien # / Accession ID		161858		nts or Addition D, Previous DSHS	

### **Laboratory Mailing Address**

Public Health Laboratory Division, MC 1947 Department of State Health Services 1100 W. 49th Street, Austin, TX 78756-3199


### **Questions About...**

**Specimen Collection/Suitability:** 512-776-7594 or 512-776-7760

Specimen Shipping: 512-776-7598 or 1-888-963-7111 ext. 7578 (toll free)

Measles Surveillance: EAIDUMeasles2025@dshs.texas.gov and Measles (Rubeola) | Texas DSHS

Submitter Accounts, Submission Forms, or Result Reports: 512-776-7578 or LabInfo@dshs.texas.gov

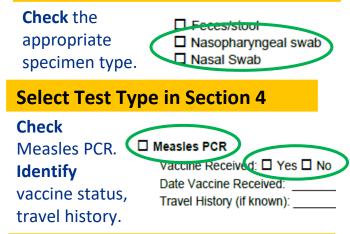

**Select Specimen Type in Section 3** 

### **Select Payor in Section 6**

**Check** the appropriate Payor. Do not leave empty.

IRED	☐ Medicaid (2) Medicaid/Medicare #:	☐ Medicare (8)
** REQUIRED	□ Submitter (3) □ BIDS (1720) □ IDEAS (1610) □ Immunizations (1609)	Private Insurance* (4)  Zoonosis (1620)  Other:

